




Hanna Instruments Inc.
Highland Industrial Park
584 Park East Drive
Woonsocket, RI 02895 USA

Technical Support for Customers
Tel. (800) 426 6287
Fax (401) 765 7575
E-mail tech@hannainst.com
www.hannainst.com


Local Sales and Customer Service Office

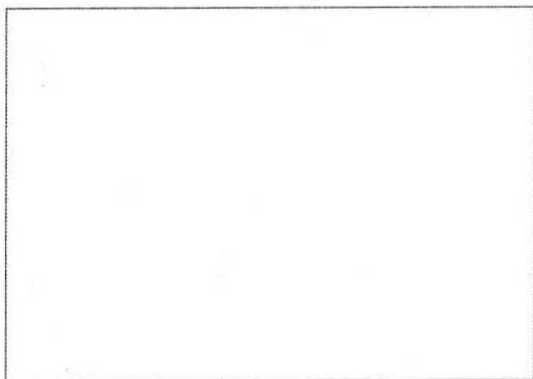

Instruction Manual

HI 2400

Dissolved Oxygen Bench Meter

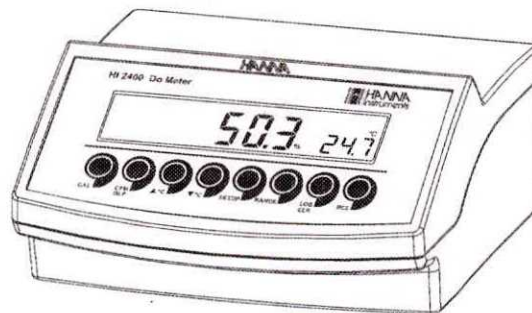

www.hannainst.com

Dear Customer,

Thank you for choosing a Hanna Instruments product.

Please read this instruction manual carefully before using these instruments.

This manual will provide you with the necessary information for correct use of these instruments, as well as a precise idea of their versatility.

If you need additional technical information, do not hesitate to e-mail us at tech@hannainst.com or view our worldwide contact list at www.hannainst.com.

WARRANTY

HI 2400 is guaranteed for two years against defects in workmanship and materials when used for its intended purpose and maintained according to instructions. Electrodes and probes are guaranteed for six months. This warranty is limited to repair or replacement free of charge.

Damage due to accidents, misuse, tampering or lack of prescribed maintenance is not covered.

If service is required, contact the dealer from whom you purchased the instrument. If under warranty, report the model number, date of purchase, serial number and the nature of the problem. If the repair is not covered by the warranty, you will be notified of the charges incurred. If the instrument is to be returned to Hanna Instruments, first obtain a Returned Goods Authorization number from the Technical Service department and then send it with shipping costs prepaid. When shipping any instrument, make sure it is properly packed for complete protection.

TABLE OF CONTENTS

WARRANTY	2
PRELIMINARY EXAMINATION	3
GENERAL DESCRIPTION	3
PROBE FUNCTIONAL DESCRIPTION	4
FUNCTIONAL DESCRIPTION	5
SPECIFICATIONS	6
OPERATIONAL GUIDE	7
DO CALIBRATION	10
GOOD LABORATORY PRACTICE (GLP)	13
SETUP	14
LOGGING	19
TEMPERATURE CALIBRATION (for technical personnel only)	22
PC INTERFACE	24
PROBE AND MEMBRANE MAINTENANCE	28
TROUBLESHOOTING GUIDE	29
ACCESSORIES	30

RECOMMENDATIONS FOR USERS

Before using this product, make sure that it is entirely suitable for the environment in which it is used.

Operation of this instrument in residential areas could cause unacceptable interferences to radio and TV equipment, requiring the operator to follow all necessary steps to correct interferences.

During operation, ESD wrist straps should be worn to avoid possible damage to the electrode by electrostatic discharges.

Any variation introduced by the user to the supplied equipment may degrade the instrument's EMC performance.

To avoid electrical shock, do not use this instrument when voltages at the measurement surface exceed 24 VAC or 60 VDC.

To avoid damage or burns, do not perform any measurement in microwave ovens.

Hanna Instruments reserves the right to modify the design, construction or appearance of its products without advance notice.

ACCESSORIES

HI 7040M	Zero Oxygen Solution, 230 ml
HI 7040L	Zero Oxygen Solution, 500 ml
HI 7041S	Refilling Electrolyte Solution, 30 ml
HI 710005	115VAC to 12VDC converter
HI 710006	230VAC to 12VDC converter
HI 76407/2	Spare probe with 2 meters (6.7') cable
HI 76407/10	Spare probe with 10 meters (33') cable
HI 76407/20	Spare probe with 20 meters (67') cable
HI 76407A/P	5 spare membranes
HI 76404N	Electrode holder

OTHER ACCESSORIES

HI 710005	Voltage adapter from 115 VAC to 12 VDC (USA plug)
HI 710006	Voltage adapter from 230 VAC to 12 VDC (European plug)
HI 710012	Voltage adapter from 240 VAC to 12 VDC (UK plug)
HI 710013	Voltage adapter from 230 VAC to 12 VDC (South African plug)
HI 710014	Voltage adapter from 230 VAC to 12 VDC (Australian plug)
ChecktempC	Pocket-size thermometer (range -50.0 to 150.0 °C)
HI 76310	Platinum 4-ring conductivity/TDS probe with temperature sensor and 1 m (3.3') cable
HI 92000	Windows® compatible software.

PRELIMINARY EXAMINATION

Remove the instrument from the packing material and examine it carefully to make sure that no damage has occurred during shipping. If there is any damage, notify your Dealer or the nearest Hanna Customer Service Center.

Each instrument is supplied with:

- HI 76407/2 DO probe with 2 m (6.7') cable
- HI 76407A membrane cap (2 pcs)
- HI 7041S electrolyte solution (30 ml)
- 12 VDC power adapter
- Instruction Manual

Note: Save all packing material until you are sure that the instrument functions correctly. All defective items must be returned in their original packing with the supplied accessories.

GENERAL DESCRIPTION

HI2400 is a logging microprocessor-based DO/Temperature bench meter. It can store up to 100 lots in memory, with up to 8000 readings. These readings can be transferred to a computer for further analysis or permanent storage.

Dissolved Oxygen is indicated in ppm (parts per million) or in %.

All measurements are automatically compensated for temperature. Salinity compensation in water allows direct determination of Dissolved Oxygen in saline waters and altitude compensation readjusts for the altitude variance.

The Dissolved Oxygen probe has a membrane covering the polarographic sensors and a built-in thermistor for temperature measurements and compensation.

This permeable PTFE membrane isolates the sensor elements from the testing solution, but allows Oxygen to pass through. When a voltage is applied across the sensor, oxygen that has passed through the membrane reacts causing a current flow, and hence determining a reading.

All rights are reserved. Reproduction in whole or in part is prohibited without the written consent of the copyright owner.

SPECIFICATIONS

RANGE	0.00 to 45.00 ppm
	0.0 to 300.0%
	0.0 to 50.0 °C (32.0 to 122°F)
RESOLUTION	0.01 ppm
	0.1%
	0.1 °C
ACCURACY @ 20 °C / 68 °F	±1.5% of full scale or ±1 digit, whichever is greater
	±0.2 °C (excluding probe error)
Typical EMC Deviation	±1.5% of full scale
	±0.5 °C
DO Calibration	Single or double point at 0% (HI 7040) and 100% (in air)
Altitude Compensation Resolution	0 to 4,000 m (13,120') 100 m (328')
Salinity Compensation Resolution	0 to 40 g/l 1 g/l
Temperature Compensation	0.0 to 50.0 °C (32.0 to 122 °F)
Probe	HI 76407/2 with 2 m (6.6') cable
Logging interval	5, 10, 30 seconds or 1, 2, 5, 10, 15, 30, 60, 120, 180 minutes
PC communication	Optoisolated USB
Power supply	12 VDC adapter
Dimensions	235 x 222 x 109 mm (9.2 x 8.7 x 4.3")
Weight	1.3 Kg (2.9 lb) kit with holder 2.1 Kg (4.6 lb)
Environment	0 – 50 °C (32 – 122 °F) max. 95% RH non condensing
Warranty	2 years

GLDxxxxff: Requests the records of the "xxx" lot number. The records are sent in frames of 10 records; "fff" is the frame number (01 first frame). (Example: Lot 13 has 53 records. The records will be sent in 6 frames, 5 with 10 records and 1 with 3 records.)

Command Parameters:

xxx - Lot number (eg: for lot number 1 xxx = 001)

fff - requested frame number - first frame is labeled 01

- Notes:**
- "Err3" is sent if the Log on demand is empty.
 - "Err4" is sent if the requested set parameter is not available.
 - "Err5" is sent if the command argument is wrong.
 - "Err6" is sent if the requested range is not available.
 - "Err7" is sent if the instrument is logging.
 - "Err8" is sent if the instrument is not in measurement mode.
 - Invalid commands will be ignored.

MDR Requests the instrument's model name and firmware code (16 ASCII chars).

GLP Requests the calibration data record.

The answer string contains:

- GLP status (1 char): represents a 4 bit hexadecimal encoding.
 - 0x04 - DO calibration available
- DO calibration data (if available), which contains:
 - the number of calibrated standards (1 char)
 - the calibration time, **yymmddhhmmss** (12 chars)
 - standards information
 - standard value, with sign and decimal point (7 chars).
 - the standard temperature, with sign and decimal point (7 chars)
 - the salinity value (4 chars)
 - the altitude value (4 chars)

PAR Requests the setup parameters setting.

The answer string contains:

- Instrument ID (4 chars)
- SETUP information (2 chars): 8 bit hexadecimal encoding.
 - 0x01 - beep ON (else OFF)
 - 0x04 - degrees Celsius (else degrees Fahrenheit)
- Log type - 01 to 12 - the corresponding interval for lot log (5 s to 180 min)
- Salinity value (4 chars)
- Altitude value (4 chars)

LLSxff: requests information about all lots and it sends the information in frames of 10 lots each (a frame contains information about 10 lots)

Command Parameters:

x = D - request DO data

f - requested frame number - first frame is labeled 01

OPERATIONAL GUIDE

POWER CONNECTION

Plug the 12 VDC adapter into the power supply socket.

- Notes:**
- This instrument use non volatile memory to retain the calibration parameters and all the other settings even when unplugged.
 - Make sure a fuse protects the main line.

PROBE CONNECTION AND PREPARATION



To take measurements, connect the DO probe to the meter securely by aligning the pins with the socket located on the back of the meter, pushing the plug in and tightening the threaded ring.

Probes shipped from Hanna Instruments are dry. To hydrate the probe and prepare it for use, connect it to the meter and proceed as follows:

1. Remove the red and black plastic cap. This cap is for shipping purposes and can be thrown away.
2. Wet the sensor by soaking the bottom 2½ cm (1") of the probe in electrolyte (HI 7041S) for 5 minutes.
3. Rinse the membrane cap (HI 76407A supplied in the kit with the meter) with electrolyte solution while shaking it gently. Refill with clean electrolyte solution.
4. Tap gently the sides of the membrane cap with your finger tip to ensure that no air bubbles are trapped. To avoid damaging the membrane, do not tap it directly on the bottom.
5. Make sure that the rubber O-ring sits properly inside the membrane cap.
6. With the sensor facing down, slowly screw the cap clockwise. Some electrolyte will overflow.

When not in use and during polarization (see page 8), use the protective transparent cap supplied in the kit with the meter.


Shipping cap

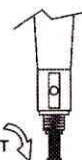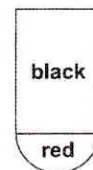

FILL FIRST

THEN TAP

THEN SCREW BACK ON

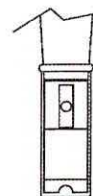

INSTRUMENT START-UP


- Turn the instrument on by pressing the ON/OFF switch.
- All LCD tags are displayed and a beep is generated while the instrument performs a self test.
- The instrument will display "Σ" blinking until initialization is complete.



- After a few seconds "Cond" message appears on the LCD to inform the user that the probe is in auto-conditioning (automatic polarization) mode.



- When this message disappears, the probe is polarized and the instrument can be calibrated.
- If the probe is disconnected, the meter will display "----".

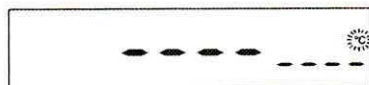

PROBE POLARIZATION

The probe is under polarization with a fixed voltage of approximately 800 mV.

Probe polarization is essential for stable measurements with the same recurring degree of accuracy.

With the probe properly polarized, oxygen is continually consumed when it passes through the sensitive diaphragm and dissolves in the electrolyte solution contained in the probe.

If polarization is interrupted, the electrolyte solution continues to be enriched with oxygen until it reaches an equilibrium with the surrounding solution.

Whenever measurements are taken with a non-polarized probe, the oxygen level revealed is both that of the tested solution, as well as that present in the electrolyte solution. This reading is incorrect.

The calibration of this instrument is very simple.

Before proceeding with the calibration, make sure the probe is ready for measurements (see page 7), i.e. the membrane cap is filled with electrolyte and the probe is connected to the meter and properly polarized.

CHR xx Change the instrument's range according to the parameter's value (xx):

- xx=06 % range
- xx=07 ppm range

The instrument will answer for these commands with:

<STX> <answer> <ETX>

where: <STX> is 02 ASCII code character (start of text)

<ETX> is 03 ASCII code character (end of text)

<answer>:

<ACK> is 06 ASCII code character (recognized command)

<NAK> is 21 ASCII code character (unrecognized command)

<CAN> is 24 ASCII code character (corrupted command)

COMMANDS REQUIRING AN ANSWER

The instrument will answer for these commands with:

<STX> <answer> <checksum> <ETX>

where the checksum is the bytes sum of the answer string sent as 2 ASCII characters.

All the answer messages are with ASCII characters.

RAS Causes the instrument to send a complete set of readings.

The answer string contains:

- Meter mode (2 chars):
 - 06 - % range
 - 07 - ppm range
- Meter status (2 chars of status byte): represents an 8 bit hexadecimal encoding.
 - 0x10 - temperature probe is connected
 - 0x01 - new GLP data available
 - 0x02 - new SETUP parameter
- Reading status (1 char): R - in range, O - over range, U - under range.
- The reading (corresponding to the selected range) - 7 ASCII chars, including sign and decimal point.
- Temperature reading - 7 ASCII chars, with sign and two decimal points, always in °C.

PC INTERFACE

Data transmission from the instrument to the PC can be done with the HI 92000 Windows[®] compatible software (optional). HI 92000 also offers graphing and an on-line help feature.

Data can be exported to the most popular spreadsheet programs for further analysis.

To connect your instrument to a PC, use a standard USB cable. Make sure that your instrument is switched off and plug one connector to the instrument's USB socket and the other to the USB port of your PC.

Note: If you are not using Hanna Instruments HI 92000 software, please see the following instructions.

SENDING COMMANDS FROM PC

It is also possible to remotely control the instrument with any terminal program. Use a standard USB cable to connect the instrument to a PC, start the terminal program and set the communication options as follows: 8, N, 1, no flow control, baud rate 9600.

COMMAND TYPES

To send a command to the instrument follow the next scheme:

< command prefix > < command > < CR >

where: < command prefix > is a selectable ASCII character between 0 and 47 (default 16).

< command > is the command code.

Note: Either small or capital letters can be used.

SIMPLE COMMANDS

RNG	Is equivalent to pressing RANGE
CAL	Is equivalent to pressing CAL
CFM	Is equivalent to pressing CFM
UPC	Is equivalent to pressing the UP arrow key
DWC	Is equivalent to pressing the DOWN arrow key
LOG	Is equivalent to pressing LOG
RCL	Is equivalent to pressing RCL
SET	Is equivalent to pressing SETUP

For an accurate calibration, it is recommended to wait at least 15 minutes to ensure precise conditioning of the probe.

Keep the protective cap on during polarization time and remove it for calibration and measurements. Follow the calibration procedure (see page 10).

SALINITY AND ALTITUDE COMPENSATION


If the sample contains significant concentration of salinity or if you are performing measurements at an altitude different from sea level, the read out values must be corrected, taking into account the lower degree of oxygen solubility in these situations (see pages 15-16).

Remember to set the altitude and/or the salinity before taking any DO measurements. The meter will automatically compensate for these factors.

DO MEASUREMENTS

Make sure that the instrument has been calibrated and the protective cap has been removed.

- Submerge the tip of the probe in the sample to be tested. Allow approximately one minute for the reading to stabilize.
- The Dissolved Oxygen value (in %) is displayed on the primary LCD and the temperature on the secondary LCD.

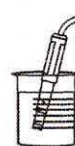

- Press RANGE to change the reading from % to ppm and vice-versa.



For accurate Dissolved Oxygen measurements, a water movement of 0.3 m/s is required. This is to ensure that the oxygen-depleted membrane surface is constantly replenished. A moving stream will provide adequate circulation.

The use of a magnetic stirrer to ensure a certain fluid velocity is recommended.

TEMPERATURE MEASUREMENTS

The probe has a built-in temperature sensor.

The measured temperature is indicated on the secondary LCD as shown above.

Allow the probe to reach thermal equilibrium before taking any measurement. This can take several minutes. The greater the difference between the temperature at which the probe was stored and the temperature of the sample, the longer the time will be.

Note: If "----" is displayed, the DO probe is not properly connected. This also indicates the possibility of a broken probe cable. If the temperature is displayed blinking, the temperature is out of range.

DO CALIBRATION

Calibrate the instrument frequently, especially if high accuracy is required.


The instrument can be calibrated in maximum 2 points: 0.0% (**zero calibration**) and 100.0% (**slope calibration**).

The zero calibration of the HI 2400 is very stable, therefore this procedure needs to be performed only whenever the probe or the membrane is replaced.


However, because the slope calibration is more critical, it is recommended to perform this procedure every week.

INITIAL PREPARATION

- Pour small quantities of HI 7040 Zero Oxygen solution into a beaker. If possible, use a plastic beaker to minimize any EMC interferences.
- Make sure the probe is ready for measurements (see probe preparation on page 7), i.e. the membrane is filled with electrolyte and the probe is connected to the meter.
- Switch the meter on by pressing the ON/OFF switch.
- For an accurate calibration, it is recommended to wait for at least 15 minutes to ensure precise conditioning of the probe.
- Remove the protective cap from the DO probe.
- Set the appropriate altitude factor (see page 15). Make sure the salinity factor is set to zero (see page 16).

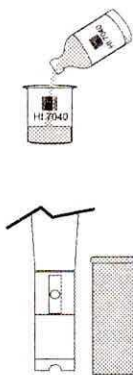

- Submerge the DO probe in the second vessel as near as possible to the reference thermometer. Allow a few seconds for the probe to stabilize.
- Use the **ARROW** keys to set the reading on the secondary LCD to that of the hot water.
- When the reading is stable, the "CFM" tag starts blinking.
- Press **CFM** to confirm. The instrument returns to measurement mode.

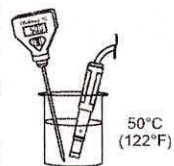

Note: If the reading is not close to the selected calibration point, "WRONG" tag will blink. In this case, check if the value set on the secondary LCD and the temperature measured with reference thermometer are close. Change the DO probe or restart calibration if necessary.

TEMPERATURE CALIBRATION (for technical personnel only)

Each meter has been factory calibrated for temperature with the supplied DO probe and is ready for measurements.

The DO probes are interchangeable and no temperature calibration is needed.


If the temperature measurements are not accurate, temperature recalibration should be performed.

For an accurate recalibration, contact your dealer or the nearest Hanna Customer Service Center, or follow the procedure below.


- Prepare a vessel containing ice and water and another one containing hot water (at a temperature of around 50 °C). Place insulation material around the vessels to minimize temperature changes.
- Use a calibrated thermometer with a resolution of 0.1 °C as a reference thermometer. Connect the HI 76407/2 DO probe to the appropriate socket.
- With the instrument off, press and hold down **CFM&SETUP**, then power on the instrument. The "CAL" tag will appear and the secondary LCD will show "0.0 °C".



- Submerge the DO probe in the vessel with ice and water as near as possible to the reference thermometer. Allow a few seconds for the probe to stabilize.

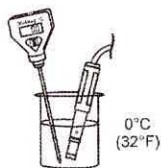


- Use the **ARROW** keys to set the reading on the secondary LCD to that of ice and water, measured by reference thermometer. When the reading is stable, the "CFM" tag starts blinking.
- Press **CFM** to confirm. The secondary LCD will show "50.0 °C".

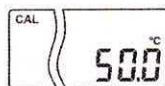

ZERO CALIBRATION

- Submerge the probe into HI 7040 zero oxygen solution and stir gently for 2-3 minutes.
- Press **CAL**. The "˘" and "⊗" tags will blink on the LCD until the reading is stable.

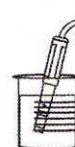

- When the reading is stable, "CFM" starts blinking. Press **CFM** to confirm the "0.0%" DO calibration.
- If the reading is within the limits ($\pm 15\%$ f.s.), the meter stores the value (and adjusts the slope point).



- Press **CAL**. The instrument will return to measurement mode and will memorize the zero calibration data. For a two-point calibration do not press **CAL** and follow the procedure below.




SLOPE CALIBRATION

It is suggested to perform the slope calibration in air.

- Rinse the probe in clean water to remove any residual zero oxygen solution.



Note: If you did not perform the zero calibration procedure, press **CAL** and then the **ARROW** keys to select the 100% DO calibration point.

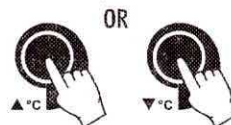

- Dry the probe tip and allow a few seconds for the LCD reading to stabilize. The "˘" and "⊗" tags will blink until the reading is stable.
- When the reading is stable, "CFM" tag starts blinking. Press **CFM** to confirm the "100.0%" DO calibration.



- If the reading is within the limits ($\pm 15\%$ f.s.), the meter stores the value (and adjusts the slope point).



- The instrument ID.

9985 In Id

SETUP

Setup mode allows viewing and modifying the following parameters:

- Salinity Factor
- Altitude Factor
- Log Interval
- Current Time (hour & minute)
- Current Date (month, day & year)
- Beep Status
- Instrument ID
- Temperature Unit

To enter the Setup mode press **SETUP** while the instrument is in measurement mode. Press **SETUP** again to exit SETUP mode.

Select a parameter with the **ARROW** keys.

Press **CAL** if you want to change a parameter value. The selected parameter will blink.

Press **RANGE** to toggle between displayed parameters.

Press the **ARROW** keys to increase or decrease the displayed value.

Press **CFM** to save the modified value or **CAL** to escape.

Press the **ARROW** keys to select the next/previous parameter.

SALINITY FACTOR

Press **CAL** when the salinity factor is displayed. The salinity factor ("0" to "40" g/l) and the "CFM" tag will start blinking.

g/L NaCl
0.5AL
SETUP

Press **ARROW** keys to change the salinity factor value.

Press **CFM** to save the modified value or press **CAL** to escape without saving.

The salinity affects the DO concentration, decreasing its value. The next table shows the maximum oxygen solubility at various temperatures and salinity levels.

LOGGING

This function allows the user to log DO (in ppm or %) together with temperature automatically, for long periods of time. All logged data can be stored into a PC through the USB port.

The memory used for storing the logged data is divided in 32 pages. The capacity of each page is 250 samples. The lot number goes from 1 to 100. The maximum capacity of the log memory is 8000 samples. Each time a new lot starts, it automatically starts from a new page. When the samples collected for a single lot reach the limit (8000 samples) or all memory pages are occupied, the meter stops logging automatically.

The appropriate logging interval can be set between 5, 10, 30 seconds or 1, 2, 5, 10, 15, 30, 60, 120, 180 minutes (see SETUP section for details).

LOGGING THE CURRENT DATA

To start the Auto LOG mode press **LOG** while in measurement mode.



When the selected interval is reached the instrument will display the current lot number on the primary LCD line, the record number on the secondary LCD line and the LOG tag (see example below: Lot 15 record 22)

LOG L 15 Log# 22

followed by the number of free records on the corresponding memory space.

LOG FREE 1228

If there are less than 6 memory locations remaining, the record number and "Lo" message will be displayed to alert the user.

LOG Lo 5

To stop the Auto LOG mode press **LOG** again. The "LOG" tag will be cleared.

Press **CFM** to save the modified value.

Press **CAL** to escape without saving.

BEEP STATUS

Press **CAL** when the beep status is displayed. The beep status ("ON" or "OFF") and "CFM" tag will start blinking.




Press the **ARROW** keys to change the beep status.

Press **CFM** to save the modified value or press **CAL** to escape without saving.

INSTRUMENT ID

Press **CAL** when "InId" is displayed. The instrument ID ("0000" to "9999") and "CFM" tag will start blinking.

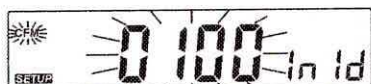

Press the **ARROW** keys to change the instrument ID value.


Press **CFM** to save the modified instrument ID value.

Press **CAL** to escape without saving.

Note: The instrument ID is downloaded to a PC as part of a logged data, set to identify its origin.

TEMPERATURE UNIT

Press **CAL** when "Unit" is displayed. The temperature unit and "CFM" tag will start blinking.

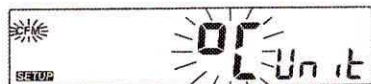

Press the **ARROW** keys to change the option.

Press **CFM** to save the modified temperature unit.

Press **CAL** to escape without saving.


°C	Salinity (g/l) at Sea Level					°F
	0 g/l	10 g/l	20 g/l	30 g/l	35 g/l	
0	14.60	13.64	12.74	11.90	11.50	32.0
2	13.81	12.91	12.07	11.29	10.91	35.6
4	13.09	12.25	11.47	10.73	10.38	39.2
6	12.44	11.65	10.91	10.22	9.89	42.8
8	11.83	11.09	10.40	9.75	9.44	46.4
10	11.28	10.58	9.93	9.32	9.03	50.0
12	10.77	10.11	9.50	8.92	8.65	53.6
14	10.29	9.68	9.10	8.55	8.30	57.2
16	9.86	9.28	8.73	8.21	7.97	60.8
18	9.45	8.90	8.39	7.90	7.66	64.4
20	9.08	8.56	8.07	7.60	7.38	68.0
22	8.73	8.23	7.77	7.33	7.12	71.6
24	8.40	7.93	7.49	7.07	6.87	75.2
25	8.24	7.79	7.36	6.95	6.75	77.0
26	8.09	7.65	7.23	6.83	6.64	78.8
28	7.81	7.38	6.98	6.61	6.42	82.4
30	7.54	7.14	6.75	6.39	6.22	86.0
32	7.29	6.90	6.54	6.19	6.03	89.6
34	7.05	6.68	6.33	6.01	5.85	93.2
36	6.82	6.47	6.14	5.83	5.68	96.8
38	6.61	6.28	5.96	5.66	5.51	100.4
40	6.41	6.09	5.79	5.50	5.36	104.0
42	6.22	5.93	5.63	5.35	5.22	107.6
44	6.04	5.77	5.48	5.21	5.09	111.2
46	5.87	5.61	5.33	5.07	4.97	114.8
48	5.70	5.47	5.20	4.95	4.85	118.4
50	5.54	5.33	5.07	4.83	4.75	122.0

Note: The relationship between salinity and chlorinity for sea water is given by the equation below:

$$\text{Salinity (g/l)} = 1.80655 \text{ Chlorinity (g/l)}$$

ALTITUDE FACTOR

Press **CAL** when the altitude factor is displayed. The altitude factor ("0" to "4000" m, in steps of 100 m; 1 meter = 3.28 feet) and the "CFM" tag will start blinking.

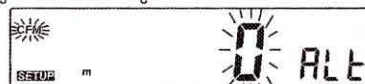


Press the **ARROW** keys to change the altitude factor value.

Press **CFM** to save the modified value or press **CAL** to escape without saving.

°C	Altitude, Meters above Sea Level															°F
	0 m	300 m	600 m	900 m	1200 m	1500 m	1800 m	2100 m	2400 m	2700 m	3000 m	3300 m	3600 m	3900 m	4000 m	
0	14.6	14.1	13.6	13.1	12.6	12.1	11.7	11.2	10.8	10.4	10.0	9.7	9.3	9.0	8.9	32.0
2	13.8	13.3	12.8	12.4	11.9	11.5	11.0	10.6	10.2	9.9	9.5	9.2	8.8	8.5	8.4	35.6
4	13.1	12.6	12.2	11.7	11.3	10.9	10.5	10.1	9.7	9.3	9.0	8.7	8.4	8.0	7.9	39.2
6	12.4	12.0	11.5	11.1	10.7	10.3	9.9	9.6	9.2	8.9	8.6	8.2	7.9	7.6	7.5	42.8
8	11.8	11.4	11.0	10.6	10.2	9.8	9.5	9.1	8.8	8.4	8.1	7.8	7.5	7.3	7.2	46.4
10	11.3	10.9	10.5	10.1	9.7	9.4	9.0	8.7	8.4	8.1	7.8	7.5	7.2	6.9	6.8	50.0
12	10.8	10.4	10.0	9.6	9.3	8.9	8.6	8.3	8.0	7.7	7.4	7.1	6.9	6.6	6.5	53.6
14	10.3	9.9	9.6	9.2	8.9	8.5	8.2	7.9	7.6	7.4	7.1	6.8	6.6	6.3	6.2	57.2
16	9.9	9.5	9.2	8.8	8.5	8.2	7.9	7.6	7.3	7.0	6.8	6.5	6.3	6.1	6.0	60.8
18	9.5	9.1	8.8	8.5	8.1	7.8	7.6	7.3	7.0	6.8	6.5	6.3	6.0	5.8	5.7	64.4
20	9.1	8.8	8.4	8.1	7.8	7.5	7.3	7.0	6.7	6.5	6.2	6.0	5.8	5.6	5.5	68.0
22	8.7	8.4	8.1	7.8	7.5	7.2	7.0	6.7	6.5	6.2	6.0	5.8	5.6	5.4	5.3	71.6
24	8.4	8.1	7.8	7.5	7.2	7.0	6.7	6.5	6.2	6.0	5.8	5.6	5.4	5.2	5.1	75.2
25	8.3	8.0	7.7	7.4	7.1	6.8	6.6	6.4	6.1	5.9	5.7	5.5	5.3	5.1	5.0	77.0
26	8.1	7.8	7.5	7.2	7.0	6.7	6.5	6.2	6.0	5.8	5.6	5.4	5.2	5.0	4.9	78.8
28	7.8	7.5	7.3	7.0	6.7	6.5	6.2	6.0	5.8	5.6	5.4	5.2	5.0	4.8	4.7	82.4
30	7.6	7.3	7.0	6.8	6.5	6.3	6.0	5.8	5.6	5.4	5.2	5.0	4.8	4.6	4.6	86.0
32	7.3	7.0	6.8	6.5	6.3	6.1	5.8	5.6	5.4	5.2	5.0	4.8	4.7	4.5	4.4	89.6
34	7.1	6.8	6.6	6.3	6.1	5.9	5.6	5.4	5.2	5.0	4.9	4.7	4.5	4.3	4.3	93.2
36	6.8	6.6	6.3	6.1	5.9	5.7	5.5	5.3	5.1	4.9	4.7	4.5	4.4	4.2	4.1	96.8
38	6.6	6.4	6.1	5.9	5.7	5.5	5.3	5.1	4.9	4.7	4.5	4.4	4.2	4.1	4.0	100.4
40	6.4	6.2	5.9	5.7	5.5	5.3	5.1	4.9	4.7	4.6	4.4	4.2	4.1	3.9	3.9	104.4
42	6.2	6.0	5.8	5.6	5.3	5.2	5.0	4.8	4.6	4.4	4.3	4.1	4.0	3.8	3.8	107.6
44	6.0	5.8	5.6	5.4	5.2	5.0	4.8	4.6	4.5	4.3	4.1	4.0	3.8	3.7	3.7	111.2
46	5.8	5.6	5.4	5.2	5.0	4.8	4.7	4.5	4.3	4.2	4.0	3.9	3.7	3.6	3.5	114.8
48	5.7	5.5	5.3	5.1	4.9	4.7	4.5	4.4	4.2	4.0	3.9	3.7	3.6	3.5	3.4	118.4
50	5.5	5.3	5.1	4.9	4.7	4.6	4.4	4.2	4.1	3.9	3.8	3.6	3.5	3.4	3.3	122.0

LOG INTERVAL

Press CAL when log interval is displayed. The log interval and "CFM" tag is displayed blinking.

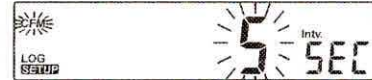

Press the ARROW keys to change the custom buffer value.

Press CFM to confirm the selection.

Press CAL to escape without saving.

CURRENT TIME

Press CAL when the current time is displayed. The hour and "CFM" tag will start blinking.



Press the ARROW keys to change the hour.

Press RANGE. The minutes will start blinking.




Press the ARROW keys to change the minutes.

Press CFM to save the modified value.

Press CAL to escape without saving.

CURRENT DATE

Press CAL when the current date is displayed. The month and "CFM" tag will start blinking.

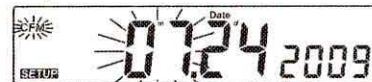

Press the ARROW keys to change the month.

Press RANGE. The day and "CFM" tag will start blinking.



Press the ARROW keys to change the day.

Press RANGE. The year and "CFM" tag will start blinking.



Press the ARROW keys to change the year.